

(12) INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(19) World Intellectual Property Organization  
International Bureau



(43) International Publication Date  
19 July 2001 (19.07.2001)

PCT

(10) International Publication Number  
**WO 01/51654 A2**

(51) International Patent Classification<sup>7</sup>: C12Q 1/68

(21) International Application Number: PCT/EP01/00306

(22) International Filing Date: 9 January 2001 (09.01.2001)

(25) Filing Language: English

(26) Publication Language: English

(30) Priority Data:  
09/481,049 11 January 2000 (11.01.2000) US

(71) Applicant: AVENTIS CROPS SCIENCE N.V. (BE/BE);  
Jozef Plateaustraat 22, B-9000 Gent (BE).

(72) Inventor: DE BEUCKELEER, Marc; Hekers 36,  
B-9052 Zwijnaarde (BE).

(81) Designated States (national): AE, AG, AL, AM, AT, AU,  
AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CR, CU, CZ,

DE, DK, DM, DZ, EE, ES, FI, GB, GD, GE, GH, GM, HR,  
HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR,  
LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ,  
NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM,  
TR, TT, TZ, UA, UG, UZ, VN, YU, ZA, ZW.

(84) Designated States (regional): ARIPO patent (GH, GM,  
KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZW), Eurasian  
patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European  
patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE,  
IT, LU, MC, NL, PT, SE, TR), OAPI patent (BF, BJ, CF,  
CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG).

**Published:**

— without international search report and to be republished  
upon receipt of that report

For two-letter codes and other abbreviations, refer to the "Guid-  
ance Notes on Codes and Abbreviations" appearing at the begin-  
ning of each regular issue of the PCT Gazette.



**WO 01/51654 A2**

(54) Title: METHODS AND KITS FOR IDENTIFYING ELITE EVENT GAT-ZM1 IN BIOLOGICAL SAMPLES

(57) Abstract: Tools are provided which allow rapid and unequivocal identification elite event GAT-ZM1 in biological samples.

Methods and kits for identifying elite event GAT-ZM1 in biological samples.

5

Background of the invention

The phenotypic expression of a transgene in a plant is determined both by the structure of the gene itself and by its location in the plant genome. At the same time the presence of the transgene (in a foreign DNA) at different locations in the genome will influence the overall phenotype of the plant in different ways. The agronomically or industrially successful introduction of a commercially interesting trait in a plant by genetic manipulation can be a lengthy procedure dependent on different factors. The actual transformation and regeneration of genetically transformed plants are only the first in a series of selection steps, which include extensive genetic characterization, breeding, and evaluation in field trials, eventually leading to the selection of an elite event.

The unequivocal identification of an elite event is becoming increasingly important in view of discussions on Novel Food/Feed, segregation of GMO and non-GMO products and the identification of proprietary material. Ideally, such identification method is both quick and simple, without the need for an extensive laboratory set-up. Furthermore, the method should provide results that allow unequivocal determination of the elite event without expert interpretation, but which hold up under expert scrutiny if necessary.

GAT-ZM1 was selected as an elite event in the development of corn resistant to the herbicide Liberty®, by transformation of corn with plasmid pUC/Ac comprising the pat gene encoding tolerance to phosphinothricin. It is commercially sold as Liberty Link® maize, such as, for instance, Liberty Link® A6460LL sold by AgriGold/Akin Seed Company. The tools for use in simple and unequivocal methods for identification elite event GAT-ZM1 in biological samples are described herein.

Summary of the invention

The present invention relates to methods for identifying elite event GAT-ZM1 in biological samples, which methods are based on primers or probes which specifically recognize the 5' and/or 3' flanking sequence of GAT-ZM1.

More specifically, the invention relates to a method comprising of amplifying a sequence of a nucleic acid present in biological samples, using a polymerase chain reaction with at least two primers, one of which recognizes the 5' or 3' flanking region of GAT-ZM1, the other which recognizes a sequence within the foreign DNA, to obtain a DNA fragment of between 100 and 350 bp. Preferably, the primers recognize a sequence within the 5' flanking region of GAT-ZM1, most preferably within the 5' flanking region of SEQ ID No. 6, and a sequence within the foreign DNA, respectively. Especially preferably, the primer recognizing the 5' flanking region comprises the nucleotide sequence of SEQ ID No. 11 and the primer recognizing a sequence within the foreign DNA comprises the nucleotide sequence of SEQ ID No. 12 described herein.

The present invention more specifically relates to a method for identifying elite event GAT-ZM1 in biological samples, which method comprises amplifying a sequence of a nucleic acid present in a biological sample, using a polymerase chain reaction with two primers having the nucleotide sequence of SEQ ID No. 11 and SEQ ID No. 12 respectively, to obtain a DNA fragment of between 180 and 220 bp, preferably of about 200 bp.

The present invention further relates to the specific flanking sequences of GAT-ZM1 described herein, which can be used to develop specific identification methods for GAT-ZM1 in biological samples. More particularly, the invention relates to the 5' and or 3' flanking regions of GAT-ZM1 which can be used for the development of specific primers and probes. The invention further relates to identification methods for the presence of GAT-ZM1 in biological samples based on the use of such specific primers or probes.

The invention further relates to kits for identifying elite event GAT-ZM1 in biological samples, said kits comprising at least one primer or probe which specifically recognizes the 5' or 3' flanking region of GAT-ZM1.

5 Preferably the kit of the invention comprises, in addition to a primer which specifically recognizes the 5' or 3' flanking region of GAT-ZM1, a second primer which specifically recognizes a sequence within the foreign DNA of GAT-ZM1, for use in a PCR identification protocol. Preferably, the kit of the invention comprises two specific primers, one of which recognizes a sequence within the 5' flanking region of GAT-ZM1, most  
10 preferably within 5' flanking region of SEQ ID No. 6, and the other which recognizes a sequence within the foreign DNA. Especially preferably, the primer recognizing the 5' flanking region comprises the nucleotide sequence of SEQ ID No. 11 and the primer recognizing the transgene comprises the nucleotide sequence of SEQ ID No. 12 described herein.

15

The invention further relates to a kit for identifying elite event GAT-ZM1 in biological samples, said kit comprising the PCR primers having the nucleotide sequence of SEQ ID No. 11 and SEQ ID No. 12 for use in the GAT-ZM1 PCR identification protocol described herein.

20

The invention also relates to a kit for identifying elite event GAT-ZM1 in biological samples, which kit comprises a specific probe having a sequence which corresponds (or is complementary to) a sequence having between 80% and 100% sequence identity with a specific region of GAT-ZM1. Preferably the sequence of the probe corresponds to a  
25 specific region comprising part of the 5' or 3' flanking region of GAT-ZM1. Most preferably the specific probe has (or is complementary to) a sequence having between 80% and 100% sequence identity to the sequence between nucleotide 286 and 466 of SEQ ID No. 6.

30 The methods and kits encompassed by the present invention can be used for different purposes such as, but not limited to the following: to identify GAT-ZM1 in plants, plant

material or in products such as, but not limited to food or feed products (fresh or processed) comprising or derived from plant material; additionally or alternatively, the methods and kits of the present invention can be used to identify transgenic plant material for purposes of segregation between transgenic and non-transgenic material; additionally  
5 or alternatively, the methods and kits of the present invention can be used to determine the quality (i.e. percentage pure material) of plant material comprising GAT-ZM1.

The invention further relates to the 5' and/or 3' flanking regions of GAT-ZM1 as well as to the specific primers and probes developed from the 5' and/or 3' flanking sequences of  
10 GAT-ZM1.

#### Detailed description

15 The incorporation of a recombinant DNA molecule in the plant genome typically results from transformation of a cell or tissue (or from another genetic manipulation). The particular site of incorporation is either due to "random" integration or is at a predetermined location (if a process of targeted integration is used).

20 The DNA introduced into the plant genome as a result of transformation of a plant cell or tissue with a recombinant DNA or "transforming DNA" is hereinafter referred to as "foreign DNA" comprising one or more "transgenes". Thus, foreign DNA may comprise both recombinant DNA as well as newly introduced, rearranged DNA of the plant. However, the term "plant DNA" in the context of the present invention will refer to DNA  
25 of the plant which is found in the same genetic locus in the corresponding wild-type plant. The foreign DNA can be characterized by the location and the configuration at the site of incorporation of the recombinant DNA molecule in the plant genome. The site in the plant genome where a recombinant DNA has been inserted is also referred to as the "insertion site" or "target site". Insertion of the recombinant DNA into the plant genome can be  
30 associated with a deletion of plant DNA, referred to as "target site deletion". A "flanking region" or "flanking sequence" as used herein refers to a sequence of at least 20 bp,

preferably at least 50 bp, and up to 5000 bp of the plant genome which is located either immediately upstream of and contiguous with or immediately downstream of and contiguous with the foreign DNA. Transformation procedures leading to random integration of the foreign DNA will result in transformants with different flanking regions, which are characteristic and unique for each transformant. When the recombinant DNA is introduced into a plant through traditional crossing, its insertion site in the plant genome, or its flanking regions will generally not be changed. An "insertion region" as used herein refers to the region corresponding to the region of at least 40 bp, preferably at least 100 bp, and up to 10000 bp, encompassed by the sequence which comprises the upstream and/or the downstream flanking region of a foreign DNA in the plant genome. Taking into consideration minor differences due to mutations within a species, an insertion region will retain, upon crossing into a plant of the same species, at least 85%, preferably 90%, more preferably 95%, and most preferably 100% sequence identity with the sequence comprising the upstream and downstream flanking regions of the foreign DNA in the plant originally obtained from transformation.

An event is defined as a (artificial) genetic locus that, as a result of genetic manipulation, carries a transgene comprising at least one copy of a gene of interest. The typical allelic states of an event are the presence or absence of the foreign DNA. An event is characterized phenotypically by the expression of the transgene. At the genetic level, an event is part of the genetic makeup of a plant. At the molecular level, an event can be characterized by the restriction map (e.g. as determined by Southern blotting), by the upstream and/or downstream flanking sequences of the transgene, the location of molecular markers and/or the molecular configuration of the transgene. Usually transformation of a plant with a transforming DNA comprising at least one gene of interest leads to a multitude of events, each of which is unique.

An elite event, as used herein, is an event which is selected from a group of events, obtained by transformation with the same transforming DNA or by back-crossing with plants obtained by such transformation, based on the expression and stability of the transgene(s) and its compatibility with optimal agronomic characteristics of the plant

comprising it. Thus the criteria for elite event selection are one or more, preferably two or more, advantageously all of the following:

- a) That the presence of the foreign DNA does not compromise other desired characteristics of the plant, such as those relating to agronomic performance or commercial value;
- 5 b) That the event is characterized by a well defined molecular configuration which is stably inherited and for which appropriate tools for identity control can be developed;
- c) That the gene(s) of interest show(s) a correct, appropriate and stable spatial and temporal phenotypic expression, both in heterozygous (or hemizygous) and homozygous condition of the event, at a commercially acceptable level in a range of environmental
- 10 conditions in which the plants carrying the event are likely to be exposed in normal agronomic use.

It is preferred that the foreign DNA is associated with a position in the plant genome that allows easy introgression into desired commercial genetic backgrounds.

15

The status of an event as an elite event is confirmed by introgression of the elite event in different relevant genetic backgrounds and observing compliance with one, two or all of the criteria e.g. a), b) and c) above.

- 20 An "elite event" thus refers to a genetic locus comprising a foreign DNA, which answers to the above-described criteria. A plant, plant material or progeny such as seeds can comprise one or more elite events in its genome.

- 25 The tools developed to identify an elite event or the plant, plant material comprising an elite event, or products which comprise plant material comprising the elite event are based on the specific genomic characteristics of the elite event, such as, a specific restriction map of the genomic region comprising the foreign DNA, molecular markers or the sequence of the flanking region(s) of the foreign DNA.

- 30 Once one or both of the flanking regions of the foreign DNA have been sequenced, primers and probes can be developed which specifically recognize this (these) sequence(s)

in the nucleic acid (DNA or RNA) of a sample by way of a molecular biological technique. For instance a PCR method can be developed to identify the elite event in biological samples (such as samples of plants, plant material or products comprising plant material). Such a PCR is based on at least two specific "primers" one recognizing a  
5 sequence within the 5' or 3' flanking region of the elite event and the other recognizing a sequence within the foreign DNA. The primers preferably have a sequence of between 15 and 35 nucleotides which under optimized PCR conditions "specifically recognize" a sequence within the 5' or 3' flanking region of the elite event and the foreign DNA of the elite event respectively, so that a specific fragment ("integration fragment") is amplified  
10 from a nucleic acid sample comprising the elite event. This means that only the targeted integration fragment, and no other sequence in the plant genome or foreign DNA, is amplified under optimized PCR conditions.

Preferably, the integration fragment has a length of between 50 and 500 nucleotides, most  
15 preferably of between 100 and 350 nucleotides. Preferably the specific primers have a sequence which is between 80 and 100% identical to a sequence within the 5' or 3' flanking region of the elite event and the foreign DNA of the elite event, respectively, provided the mismatches still allow specific identification of the elite event with these primers under optimized PCR conditions. The range of allowable mismatches however,  
20 can easily be determined experimentally and are known to a person skilled in the art.

As the sequence of the primers and their relative location in the genome are unique for the elite event, amplification of the integration fragment will occur only in biological samples comprising (the nucleic acid of) the elite event. Preferably when performing a PCR to  
25 identify the presence of GAT-ZM1 in unknown samples, a control is included of a set of primers with which a fragment within a "housekeeping gene" of the plant species of the event can be amplified. Housekeeping genes are genes that are expressed in most cell types and which are concerned with basic metabolic activities common to all cells. Preferably, the fragment amplified from the housekeeping gene is a fragment which is  
30 larger than the amplified integration fragment. Depending on the samples to be analyzed, other controls can be included.



Standard PCR protocols are described in the art, such as in "PCR Applications Manual" (Roche Molecular Biochemicals, 2nd Edition, 1999). The optimal conditions for the PCR, including the sequence of the specific primers, is specified in a "PCR identification protocol" for each elite event. It is however understood that a number of parameters in the PCR identification protocol may need to be adjusted to specific laboratory conditions, and may be modified slightly to obtain similar results. For instance, use of a different method for preparation of DNA may require adjustment of, for instance, the amount of primers, polymerase and annealing conditions used. Similarly, the selection of other primers may dictate other optimal conditions for the PCR identification protocol. These adjustments will however be apparent to a person skilled in the art, and are furthermore detailed in current PCR application manuals such as the one cited above.

Alternatively, specific primers can be used to amplify an integration fragment that can be used as a "specific probe" for identifying GAT-ZM1 in biological samples. Contacting nucleic acid of a biological sample, with the probe, under conditions which allow hybridization of the probe with its corresponding fragment in the nucleic acid, results in the formation of a nucleic acid/probe hybrid. The formation of this hybrid can be detected (e.g. labeling of the nucleic acid or probe), whereby the formation of this hybrid indicates the presence of GAT-ZM1. Such identification methods based on hybridization with a specific probe (either on a solid phase carrier or in solution) have been described in the art. The specific probe is preferably a sequence which, under optimized conditions, hybridizes specifically to a region within the 5' or 3' flanking region of the elite event and preferably also comprising part of the foreign DNA contiguous therewith (hereinafter referred to as "specific region"). Preferably, the specific probe comprises a sequence of between 50 and 500 bp, preferably of 100 to 350 bp which is at least 80%, preferably between 80 and 85%, more preferably between 85 and 90%, especially preferably between 90 and 95%, most preferably between 95% and 100% identical (or complementary) to the nucleotide sequence of a specific region. Preferably, the specific probe will comprise a sequence of about 15 to about 100 contiguous nucleotides identical (or complementary) to a specific region of the elite event.

A "kit" as used herein refers to a set of reagents for the purpose of performing the method of the invention, more particularly, the identification of the elite event GAT-ZM1 in biological samples. More particularly, a preferred embodiment of the kit of the invention comprises at least one or two specific primers, as described above. Optionally, the kit can  
5 further comprise any other reagent described herein in the PCR identification protocol. Alternatively, according to another embodiment of this invention, the kit can comprise a specific probe, as described above, which specifically hybridizes with nucleic acid of biological samples to identify the presence of GAT-ZM1 therein. Optionally, the kit can further comprise any other reagent (such as but not limited to hybridizing buffer, label) for  
10 identification of GAT-ZM1 in biological samples, using the specific probe.

The kit of the invention can be used, and its components can be specifically adjusted, for purposes of quality control (e.g., purity of seed lots), detection of the elite event in plant material or material comprising or derived from plant material, such as but not limited to  
15 food or feed products.

As used herein, "sequence identity" with regard to nucleotide sequences (DNA or RNA), refers to the number of positions with identical nucleotides divided by the number of nucleotides in the shorter of the two sequences. The alignment of the two nucleotide  
20 sequences is performed by the Wilbur and Lipmann algorithm (Wilbur and Lipmann, 1983, Proc. Nat. Acad. Sci. USA 80:726) using a window-size of 20 nucleotides, a word length of 4 nucleotides, and a gap penalty of 4. Computer-assisted analysis and interpretation of sequence data, including sequence alignment as described above, can, e.g., be conveniently performed using the programs of the Intelligenetics<sup>TM</sup> Suite  
25 (Intelligenetics Inc., CA) or the sequence analysis software package of the Genetics Computer Group (GCG, University of Wisconsin Biotechnology center). Sequences are indicated as "essentially similar" when such sequences have a sequence identity of at least about 75%, particularly at least about 80%, more particularly at least about 85%, quite particularly about 90%, especially about 95%, more especially about 100%. It is clear than  
30 when RNA sequences are said to be essentially similar or have a certain degree of sequence identity with DNA sequences, thymidine (T) in the DNA sequence is considered

equal to uracil (U) in the RNA sequence. "Complementary to" as used herein refers to the complementarity between the A and T (U), and G and C nucleotides in nucleotide sequences.

- 5 The term "primer" as used herein encompasses any nucleic acid that is capable of priming the synthesis of a nascent nucleic acid in a template-dependent process, such as PCR. Typically, primers are oligonucleotides from 10 to 30 basepairs, but longer sequences can be employed. Primers may be provided in double-stranded form, though the single-stranded form is preferred. Probes can be used as primers, but are designed to bind to the target DNA or RNA and need not be used in an amplification process.
- 10

- The term "recognizing" as used herein when referring to specific primers, refers to the fact that the specific primers specifically hybridize to a nucleic acid sequence in the elite event under the conditions set forth in the method (such as the conditions of the PCR identification protocol), whereby the specificity is determined by the presence of positive and negative controls.
- 15

- The term "hybridizing" as used herein when referring to specific probes, refer to the fact that the probe binds to a specific region in the nucleic acid sequence of the elite event under standard stringency conditions. Standard stringency conditions as used herein refers to the conditions for hybridization described herein or to the conventional hybridizing conditions as described by Sambrook et al., 1989 (Molecular Cloning: A Laboratory Manual, Second Edition, Cold Spring Harbour Laboratory Press, NY) which for instance can comprise the following steps: 1) immobilizing plant genomic DNA fragments on a filter, 2) prehybridizing the filter for 1 to 2 hours at 42°C in 50% formamide, 5 X SSPE, 2 X Denhardt's reagent and 0.1% SDS, or for 1 to 2 hours at 68°C in 6 X SSC, 2 X Denhardt's reagent and 0.1% SDS, 3) adding the hybridization probe which has been labeled, 4) incubating for 16 to 24 hours, 5) washing the filter for 20 min. at room temperature in 1X SSC, 0.1 %SDS, 6) washing the filter three times for 20 min. each at 68°C in 0.2 X SSC, 0.1 %SDS, and 7) exposing the filter for 24 to 48 hours to X-ray film at -70°C with an intensifying screen.
- 20
- 25
- 30

As used in herein, a biological samples is a sample of a plant, plant material or products comprising plant material. The term "plant" is intended to encompass corn (*Zea mays*) plant tissues, at any stage of maturity, as well as any cells, tissues, or organs taken from or  
5 derived from any such plant, including without limitation, any seeds, leaves, stems, flowers, roots, single cells, gametes, cell cultures, tissue cultures or protoplasts. "Plant material", as used herein refers to material which is obtained or derived from a plant. Products comprising plant material relate to food, feed or other products which are produced using plant material or can be contaminated by plant material. It is understood  
10 that, in the context of the present invention, such biological samples are tested for the presence of nucleic acids specific for GAT-ZM1, implying the presence of nucleic acids in the samples. Thus the methods referred to herein for identifying elite event GAT-ZM1 in biological samples, relate to the identification in biological samples of nucleic acids which comprise the elite event.

15 As used herein "comprising" is to be interpreted as specifying the presence of the stated features, integers, steps, reagents or components as referred to, but does not preclude the presence or addition of one or more features, integers, steps or components, or groups thereof. Thus, e.g., a nucleic acid or protein comprising a sequence of nucleotides or  
20 amino acids, may comprise more nucleotides or amino acids than the actually cited ones, i.e., be embedded in a larger nucleic acid or protein. A chimeric gene comprising a DNA sequence which is functionally or structurally defined, may comprise additional DNA sequences, etc.

25 The following examples describe the identification of the development of tools for the identification of elite event GAT-ZM1 in biological samples.

Unless otherwise stated, all recombinant DNA techniques are carried out according to standard protocols as described in Sambrook et al. (1989) *Molecular Cloning: A*  
30 *Laboratory Manual*, Second Edition, Cold Spring Harbour Laboratory Press, NY and in Volumes 1 and 2 of Ausubel et al. (1994) *Current Protocols in Molecular Biology*,

Current Protocols, USA. Standard materials and methods for plant molecular work are described in Plant Molecular Biology Labfax (1993) by R.D.D. Croy published by BIOS Scientific Publications Ltd (UK) and Blackwell Scientific Publications, UK.

5

In the description and examples, reference is made to the following sequences:

- |                   |   |
|-------------------|---|
| SEQ ID No. 1:     | sequence of the genetic elements of vector pUC/Ac                 |
| SEQ ID No. 2:     | primer MDB286   |
| 10 SEQ ID No. 3:  | primer MDB391   |
| SEQ ID No. 4:     | primer MDB411   |
| SEQ ID No. 5:     | primer MDB420   |
| SEQ ID No. 6:     | nucleotide sequence comprising a 5' flanking region of<br>GAT-ZM1 |
| 15 SEQ ID No. 7:  | primer MDB439   |
| SEQ ID No. 8:     | primer VDS44  |
| SEQ ID No. 9:     | primer MDB522   |
| SEQ ID No. 10:    | nucleotide sequence comprising a 3' flanking region of<br>GAT-ZM1 |
| 20 SEQ ID No. 11: | primer COR17  |
| SEQ ID No. 12:    | primer COR18  |
| SEQ ID No. 13:    | primer COR15  |
| SEQ ID No. 14:    | primer COR16  |

25

Brief description of the drawing

The following Examples, not intended to limit the invention to specific embodiments described, may be understood in conjunction with the accompanying Figure, incorporated  
 5 herein by reference, in which:

**Fig. 1: Scoring of unknowns using the PCR Identification protocol developed for GAT-ZM1.** Loading sequence of the gel: Unknowns: lanes 1, 2, 5, 6, 8, 11, 13, 14, DNA samples from corn plants comprising the transgenic event GAT-ZM1; lanes  
 10 4, 9, 10, 12, DNA samples from a corn plants not comprising elite event GAT-ZM1; lane 3, PCR failure. Control lanes: lanes 19, 21, control DNA samples from corn plants comprising elite event GAT-ZM1; lanes 20, 22, control DNA samples from wild-type corn plants; lane 23, no template control; lane 24, molecular weight marker.

15

Examples**1. Identification of the flanking regions of elite event GAT-ZM1**

20 Herbicide-resistant corn was developed by transformation of corn with the pUC/Ac vector comprising the coding sequence of a *pat* gene encoding the enzyme phosphinothricin-acetyl-transferase, under the control of the constitutive 35S promoter from Cauliflower Mosaic virus. A detailed description of the genetic elements of pUC/Ac is provided in Table 1. The nucleotide sequence of the genetic elements of pUC/Ac is provided in SEQ  
 25 ID No. 1.

Table: genetic elements of vector pUC/Ac

Nucleotide no.	Genetic element
412 – 618	<b>35S terminator from Cauliflower Mosaic Virus</b> from the vector pDH51(Pietrzak M. et al., Nucl. Acids Res. 14, (1986), pp. 5857-5868)

619 – 636	Synthetic polylinker sequences
637 – 1188	<b>Synthetic pat gene</b> (amino acid sequence from <i>Streptomyces viridochromogenes</i> ) (Strauch et al. (1993) European patent 275957 B1)
1189 – 1216	Synthetic polylinker sequences
1217 – 1746	<b>35S promoter from Cauliflower Mosaic Virus</b> from the vector pDH51(Pietrzak M. et al., 1986)
1747 – 411	Sequence of the vector pUC18, including the $\beta$ -lactamase gene (pos. 2923 – 3783) and the origin of replication at pos. 2164. (Yanisch-Perron et al., (1985), Gene 33, pp. 103-119)

Elite event GAT-ZM1 was selected based on an extensive selection procedure based on good expression and stability of the herbicide resistance gene and its compatibility with optimal agronomic characteristics.

The sequence of the regions flanking the foreign DNA in the GAT-ZM1 event was determined using the thermal asymmetric interlaced (TAIL-) PCR method described by Liu et al. (1995, Plant J. 8(3):457-463). This method utilizes three nested primers in successive reactions together with a shorter arbitrary degenerate primer so that the relative amplification efficiencies of specific and non-specific products can be thermally controlled. The specific primers were selected for annealing to the border of the foreign DNA and based on their annealing conditions. A small amount (5 $\mu$ l) of unpurified secondary and tertiary PCR products were analyzed on a 1% agarose gel. The tertiary PCR product was used for preparative amplification, purified and sequenced on an automated sequencer using the DyeDeoxy Terminator cycle kit.

#### 1.1. Right (5') flanking region

The primers used were:

	Sequence (5' → 3')	Position in pUC/Ac
Degenerate primer MDB286	NgT.CgA.SWg.ANA.WgA.A (SEQ ID No. 2)	-
Primary TAIL MDB391	Tgg.ATA.CAA.gCA.Tgg.Tgg.ATg.g (SEQ ID No. 3)	715←736
Secondary TAIL MDB411	Agg.CAT.gCC.gCT.gAA.ATC.ACC (SEQ ID NO. 4)	606←626
Tertiary TAIL MDB420	GgT.TTC.gCT.CAT.gTg.TTg.AgC (SEQ ID No. 5)	507←527

Whereby: N = A,C,T or g; S = C or g; W = A or T

The fragment amplified using MDB286-MDB420 was ca. 1100 bp, the complete sequence of which was determined (SEQ ID No. 6). The sequence between nucleotide 1 and 341 corresponds to plant DNA, while the sequence between nucleotide 342 and 1041 corresponds to T-DNA.

## 1.2. Left (3') flanking region

10 The primers used were:

	Sequence (5' → 3')	Position in pUC/Ac
Degenerate primer MDB286	NgT.CgA.SWg.ANA.WgA.A (SEQ ID No. 2)	-
Primary TAIL MDB439	CTC.ATg.gTT.ATg.gCA.gCA.CTg.C (SEQ ID No. 7)	3401-3422
Secondary TAIL VDS44	CTg.TCA.TgC.CAT.CCg.TAA.gAT.gC (SEQ ID No. 8)	3435-3457
Tertiary TAIL MDB522	TgC.TTT.gAA.gAC.gTg.gTT.gg (SEQ ID No. 9)	1326-1345

Whereby: N = A,C,T or g; S = C or g; W = A or T

The fragment amplified using MDB286-MDB522 was ca. 450 bp, the complete sequence of which was determined (SEQ ID No. 10). The sequence between nucleotide 1 and 342



corresponds to T-DNA, while the sequence between nucleotide 343 and 484 corresponds to plant DNA.

## 2. Development of a Polymerase Chain reaction identification protocol

### 2.1. Primers

Specific primers were developed which recognize sequences within the elite event. More particularly, a primer was developed which recognizes a sequence within the 5' flanking region of GAT-ZM1. A second primer was then selected within the sequence of the foreign DNA so that the primers span a sequence of about 200 bp. The following primers were found to give particularly clear and reproducible results in a PCR reaction on GAT-ZM1 DNA:

COR17: 5'-ggg.TgA.gCT.CgA.ATg.TTg.TTC.T-3' (SEQ ID 11)  
(target: plant DNA)

COR18: 5'-TCT.TAg.ACg.TCA.ggT.ggC.ACT.T-3' (SEQ ID 12)  
(target: T-DNA)

Primers targeting an endogenous sequence are preferably included in the PCR cocktail. These primers serve as an internal control in unknown samples and in the DNA positive control. A positive result with the endogenous primer-pair demonstrates that there is ample DNA of adequate quality in the genomic DNA preparation for a PCR product to be generated. The endogenous primers were selected to recognize a housekeeping gene in *Zea mays*:

COR15: 5'-AgC.gTC.AAg.gAT.CAT.Tgg.TgT.C-3' (SEQ ID 13)  
(located in *Zea Mays* alcohol dehydrogenase 1 gene (X04050))

COR16: 5'-ggC.CAA.gTT.CAg.CAT.AAg.CTg.T-3' (SEQ ID 14)  
(located in *Zea Mays* alcohol dehydrogenase 1 gene (X04050))

## 5 2.2. Amplified fragments

The expected amplified fragments in the PCR reaction are:

	For primer pair COR15-COR16:	513bp (endogenous control)
10	For primer pair COR17-COR18:	202bp (GAT-ZM1 elite Event)

## 2.3. Template DNA

15 Template DNA was prepared from a leaf punch according to Edwards *et al.* (Nucleic Acid Research, 19, p1349, 1991). When using DNA prepared with other methods, a test run utilizing different amounts of template should be done. Usually 50 ng of genomic template DNA yields the best results.

20

## 2.4. Assigned positive and negative controls

To avoid false positives or negatives, it was determined that the following positive and negative controls should be included in a PCR run:

25

- Master Mix control (DNA negative control). This is a PCR in which no DNA is added to the reaction. When the expected result, no PCR products, is observed this indicates that the PCR cocktail was not contaminated with target DNA.

- A DNA positive control (genomic DNA sample known to contain the transgenic sequences). Successful amplification of this positive control demonstrates that the PCR was run under conditions which allow for the amplification of target sequences.
- 5 - A wildtype DNA control. This is a PCR in which the template DNA provided is genomic DNA prepared from a non-transgenic plant. When the expected result, no amplification of a transgene PCR product but amplification of the endogenous PCR product, is observed this indicates that there is no detectable transgene background amplification in a genomic DNA sample.

10

## 2.5. PCR conditions

Optimal results were obtained under the following conditions:

15

- the PCR mix for 25 $\mu$ l reactions contains:

2.5  $\mu$ l template DNA  
2.5  $\mu$ l 10x Amplification Buffer (supplied with Taq polymerase)  
20 0.5  $\mu$ l 10 mM dNTP's  
0.5  $\mu$ l COR17 (10pmoles/ $\mu$ l)  
0.5  $\mu$ l COR18 (10pmoles/ $\mu$ l)  
0.25  $\mu$ l COR15 (10pmoles/ $\mu$ l)  
0.25  $\mu$ l COR16 (10pmoles/ $\mu$ l)  
25 0.1  $\mu$ l Taq DNA polymerase (5 units/ $\mu$ l)  
water up to 25  $\mu$ l

- the thermocycling profile to be followed for optimal results is the following:

30

4 min. at 95°C

Followed by: 1 min. at 95°C  
1 min. at 57°C  
2 min. at 72°C  
For 5 cycles

5

Followed by: 30 sec. at 92°C  
30 sec. at 57°C  
1 min. at 72°C  
For 25 cycles

10

Followed by: 5 minutes at 72°C

## 15 2.6. Agarose gel analysis

To optimally visualise the results of the PCR it was determined that between 10 and 20µl of the PCR samples should be applied on a 1.5% agarose gel (Tris-borate buffer) with an appropriate molecular weight marker (e.g. 100bp ladder PHARMACIA).

20

## 2.7. Validation of the results

It was determined that data from transgenic plant DNA samples within a single PCR run and a single PCR cocktail should not be acceptable unless 1) the DNA positive control shows the expected PCR products (transgenic and endogenous fragments), 2) the DNA negative control is negative for PCR amplification (no fragments) and 3) the wild-type DNA control shows the expected result (endogenous fragment amplification).

30 When following the PCR Identification Protocol for GAT-ZM1 as described above, lanes showing visible amounts of the transgenic and endogenous PCR products of the expected

sizes, indicate that the corresponding plant from which the genomic template DNA was prepared, has inherited the GAT-ZM1 elite event. Lanes not showing visible amounts of either of the transgenic PCR products and showing visible amounts of the endogenous PCR product, indicate that the corresponding plant from which the genomic template DNA was prepared, does not comprise the elite event. Lanes not showing visible amounts of the endogenous and transgenic PCR products, indicate that the quality and/or quantity of the genomic DNA didn't allow for a PCR product to be generated. These plants cannot be scored. The genomic DNA preparation should be repeated and a new PCR run, with the appropriate controls, has to be performed.

10

## 2.8. Use of discriminating PCR protocol to identify GAT-ZM1

Before attempting to screen unknowns, a test run, with all appropriate controls, has to be performed. The developed protocol might require optimization for components that may differ between labs (template DNA preparation, Taq DNA polymerase, quality of the primers, dNTP's, thermocycler, etc.).

Amplification of the endogenous sequence plays a key role in the protocol. One has to attain PCR and thermocycling conditions that amplify equimolar quantities of both the endogenous and transgenic sequence in a known transgenic genomic DNA template. Whenever the targeted endogenous fragment is not amplified or whenever the targeted sequences are not amplified with the same ethidium bromide staining intensities, as judged by agarose gel electrophoresis, optimization of the PCR conditions may be required.

*Zea mays* leaf material from a number of plants, some of which comprising GAT-ZM1 were tested according to the above-described protocol. Samples from elite event GAT-ZM1 and from *Zea mays* wild-type were taken as positive and negative controls, respectively.

Figure 1 illustrates the result obtained with the elite event PCR identification protocol for GAT-ZM1 on a number of corn plant samples (lanes 1 to 14). The samples in lanes 1, 2, 5,

6, 7, 8, 11, 13, and 14 were found to contain the elite event as the 202 bp band is detected, while the samples in lanes 4, 9, 10, and 12 do not comprise GAT-ZM1. Lane 3 indicates a PCR failure, as the control band is not detected. Lanes 19 and 20 represent GAT-ZM1 positive control samples, lanes 20 and 22 represent non-transgenic *Zea mays* controls; lane 5 23 represents the negative control (water) sample, and lane 24 the Molecular Weight Marker (100 bp).

### 3. Use of a specific integration fragment as a probe for detection of material comprising GAT-ZM1

10

A specific integration fragment of GAT-ZM1 is obtained by PCR amplification using specific primers COR17 (SEQ ID No. 11) and COR18 (SEQ ID No. 12) or by chemical synthesis and is labeled. This integration fragment is used as a specific probe for the detection of GAT-ZM1 in biological samples. Nucleic acid is extracted from the samples 15 according to standard procedures. This nucleic acid is then contacted with the specific probe under hybridization conditions which are optimized to allow formation of a hybrid. The formation of the hybrid is then detected to indicate the presence of GAT-ZM1 nucleic acid in the sample. Optionally, the nucleic acid in the samples is amplified using the specific primers prior to contact with the specific probe. Alternatively, the nucleic acid is 20 labeled prior to contact with the specific probe instead of the integration fragment. Optionally, the specific probe is attached to a solid carrier (such as, but not limited to a filter, strip or beads), prior to contact with the samples.

Claims

1. A method for identifying elite event GAT-ZM1 in biological samples, which method comprises detection of a GAT-ZM1 specific region with a specific primer or probe  
5 which specifically recognizes the 5' or 3' flanking region of GAT-ZM1.
2. The method of claim 1, said method comprising amplifying a DNA fragment of between 100 and 350 bp from a nucleic acid present in said biological samples using a polymerase chain reaction with at least two primers, one of which recognizes the  
10 5' or 3' flanking region of GAT-ZM1, the other which recognizes a sequence within the foreign DNA.
3. The method of claim 2, wherein one of said primers recognizes a sequence within the 5' flanking region of GAT-ZM1, and said other primer recognizes a sequence  
15 within the foreign DNA.
4. The method of claim 3, wherein one of said primers recognizes a sequence within the 5' flanking region of SEQ ID No. 6, and the other recognizes a sequence within the foreign DNA.  
20
5. The method of claim 4, wherein said primers comprise the sequence of SEQ ID No. 11 and SEQ ID No. 12, respectively.
6. The method of claim 2, which method comprises amplifying a fragment of between  
25 150 and 220 bp using the GAT-ZM1 identification protocol, whereby the sequence of said primers corresponds to the nucleotide sequence of SEQ ID No 11 and SEQ ID No 12 respectively.
7. The method of claim 6, which method comprises amplifying a fragment of about  
30 202 bp, using the GAT-ZM1 identification protocol.

8. A kit for identifying elite event GAT-ZM1 in biological samples, said kit comprising at least one PCR primer, which recognizes a sequence within the 3' or 5' border flanking region of GAT-ZM1.
- 5 9. The kit of Claim 8, which further comprises at least a second PCR primer which recognizes a sequence within the foreign DNA of GAT-ZM1.
- 10 10. The kit of claim 8, wherein said at least one PCR primer recognizes a sequence within the 5' flanking region of SEQ ID No. 6.
- 10 11. The kit of claim 8, wherein said at least two PCR primers comprise the sequence of SEQ ID No. 11 and SEQ ID No. 12, respectively.
- 15 12. A primer for use in a GAT-ZM1 PCR identification protocol, having a sequence which, under optimized PCR conditions specifically recognizes a sequence within the 5' or 3' flanking region of GAT-ZM1.
- 20 13. The primer of claim 12, having a sequence which has at least 80% sequence identity with a sequence within SEQ ID No. 6 or SEQ ID No. 10.
- 20 14. The primer having the sequence of SEQ ID No. 11.
- 25 15. The primer having the sequence of SEQ ID No. 12.
- 25 16. The method of claim 1, which method comprising hybridizing a nucleic acid of biological samples with a specific probe for GAT-ZM1.
- 30 17. The method of claim 16, wherein the sequence of said specific probe has at least 80% sequence identity with a sequence comprising part of the 5' flanking sequence of GAT-ZM1 and the sequence of the foreign DNA contiguous therewith.



18. The method of claim 17, wherein the sequence of said specific probe has at least 80% sequence identity with SEQ ID No. 6, from nucleotide 286 to 487.
19. A kit for identifying elite event GAT-ZM1 in biological samples, said kit comprising  
5 a specific probe, capable of hybridizing specifically to a specific region of GAT-ZM1.
20. The kit of claim 19, wherein the sequence of said specific probe has at least 80% sequence identity with a sequence comprising part of the 5' flanking sequence of  
10 GAT-ZM1 and the sequence of the foreign DNA contiguous therewith.
21. The kit of claim 20, wherein the sequence of said specific probe has at least 80% sequence identity with SEQ ID No. 6, from nucleotide 286 to 487, or the complement thereof.
- 15 22. A specific probe for the identification of elite event GAT-ZM1 in biological samples.
23. The probe of claim 22, which has at least 80% sequence identity with a sequence  
20 comprising part of the 5' flanking sequence of GAT-ZM1 and the sequence of the foreign DNA contiguous therewith, or the complement thereof.
24. The probe of claim 23 which has at least 80% sequence identity with SEQ ID No. 6, from nucleotide 286 to 487 or the complement thereof.
- 25 25. A specific probe probe for the identification of elite event GAT-ZM1 in biological samples, the sequence of being essentially similar to SEQ ID No. 6, from nucleotide 286 to 487 or the complement thereof.

26. A method for confirming seed purity, which method comprises detection of a GAT-ZM1 specific region with a specific primer or probe which specifically recognizes the 5' or 3' flanking region of GAT-ZM1, in seed samples.
- 5 27. A method for screening seeds for the presence of GAT-ZM1, which method comprises detection of a GAT-ZM1 specific region with a specific primer or probe which specifically recognizes the 5' or 3' flanking region of GAT-ZM1, in samples of seed lots.

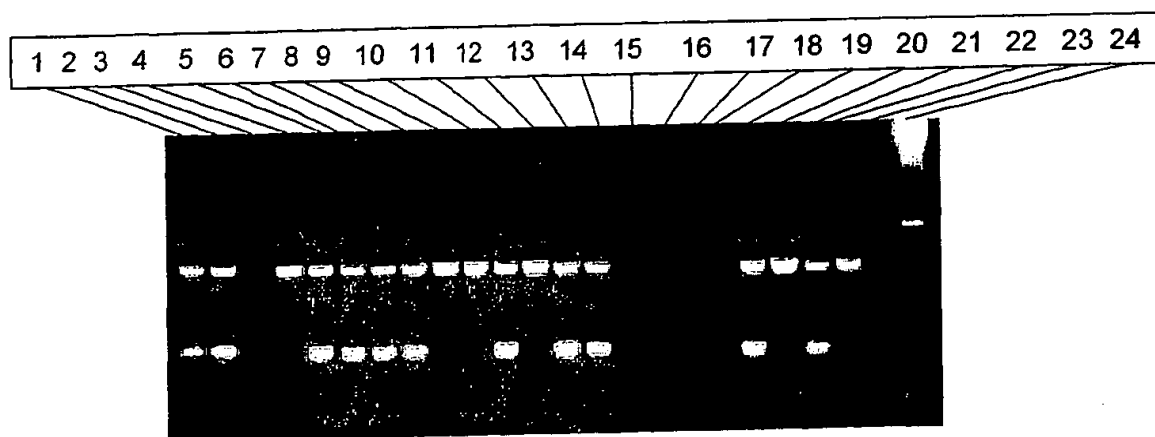


Figure 1

## SEQUENCE LISTING

<110> AVENTIS CROPSCIENCE N.V.

<120> Methods and kits for identifying elite event GAT-ZM1  
in biological samples

<130> EE-ZM1

<140>

<141>

<150> US 09/481049

<151> 2000-01-11

<160> 14

<170> PatentIn Ver. 2.0

<210> 1

<211> 3983

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: genetic  
elements of pUC/Ac

<220>

<221> terminator

<222> (412)..(618)

<223> 35S terminator from Cauliflower Mosaic Virus

<220>

<221> misc\_feature

<222> (619)..(636)

<223> synthetic polylinker sequence

<220>

<221> gene

<222> (637)..(1188)

<223> Synthetic pat-gene coding sequence

<220>

<221> promoter

<222> (1217)..(1746)

<223> 35S promoter from Cauliflower Mosaic Virus

&lt;400&gt; 1

tcgcgcggttt cggatgatgac ggtgaaaacc tctgacacat gcagctcccg gagacgggtca 60  
cagcttgtct gtaagcggat gccgggagca gacaagcccg tcagggcgcg tcagcgggtg 120  
ttggcgggtg tcggggctgg cttaactatg cggcatcaga gcagattgta ctgagagtgc 180  
accatattgcg gtgtgaaata ccgcacagat gcgtaaggag aaaataccgc atcaggcgcc 240  
attcgccatt caggctgctc aactgttggg aaggcgatc ggtgcgggccc tcttcgctat 300  
tacgccagct ggcgaaaggg ggatgtgctg caaggcgatt aagttgggta acgccagggt 360  
tttcccagtc acgacgttgt aaaacgacgg ccagtgccaa gcttgaattc gagctcggta 420  
cccactggat tttggtttta ggaattagaa attttattga tagaagtatt ttacaaatac 480  
aaatacatac taagggtttc ttatatgctc aacacatgag cgaaacccta taagaaccct 540  
aattccctta tctgggaact actcacacat tattatagag agagatagat ttgtagagag 600  
agactggtga tttcagcggc atgcctgcag gtcgactcag atctgggtaa ctggcctaac 660  
tggccttga ggagctggca actcaaaatc cctttgccaa aaaccaacat catgccatcc 720  
accatgcttg tatccagctg cgcgcaatgt acccggggct gtgtatccca aagcctcatg 780  
caacctaaca gatggatcgt ttggaaggcc tataacagca accacagact taaaaccttg 840  
cgcctccata gacttaagca aatgtgtgta caatgtggat cctaggccca acctttgatg 900  
cctatgtgac acgtaaacag tactctcaac tgtccaatcg taagcgttcc tagccttcca 960  
gggcccagcg taagcaatac cagccacaac accctcaacc tcagcaacca accaagggtg 1020  
tctatcttgc aacctctcta gatcatcaat ccactcttgt ggtgtttgtg gctctgtcct 1080  
aaagttcact gtagacgtct caatgtaatg gttaacgata tcacaaaccg cggccatatt 1140  
agctgctgta gctggcctaa tctcaactgg tctcctctcc ggagacatgt cgactctaga 1200  
ggatccccgg gtaccctgtc ctctccaaat gaaatgaact tccttatata gaggaagggt 1260  
cttgcgagg atagtgggat tgtgcgtcat cccttacgtc agtggagata tcacatcaat 1320  
ccacttgctt tgaagacgtg gttggaacgt cttctttttc cacgatgctc ctctgtgggtg 1380  
ggggtccatc tttgggacca ctgtcggcag aggcattctc aacgatggcc tttcctttat 1440  
cgcaatgatg gcatttgtag gagccacctt ccttttccac tatcttcaca ataaagtgc 1500  
agatagctgg gcaatggaat ccgaggaggt ttccggatat taccctttgt tgaaaagtct 1560  
caattgccct ttggtcttct gagactgtat ctttgatatt tttggagtag acaagcgtgt 1620  
cgtgctccac catgttgacg aagattttct tctgtcatt gagtcgtaag agactctgta 1680  
tgaactgttc gccagtcttt acggcgagtt ctgttaggtc ctctatttga atctttgact 1740  
ccatgggaat tcgtaatcat ggtcatagct gtttctctgt tgaaattgtt atccgctcac 1800  
aattccacac aacatacgag ccggaagcat aaagtgtaaa gcctgggggtg cctaattgagt 1860  
gagctaactc acattaattg cgttgcgctc actgcccgct ttccagtcgg gaaacctgtc 1920  
gtgccagctg cattaatgaa tcggccaacg cgcggggaga ggcggtttgc gtattgggag 1980  
ctcttccgct tcctcgctca ctgactcgct gcgctcggtc gttcggctgc ggcgagcggg 2040  
atcagctcac tcaaaggcgg taatacgggt atccacagaa tcaggggata acgcaggaaa 2100  
gaacatgtga gcaaaaggcc agcaaaaggc caggaaaccgt aaaaaggccg cgttgctggc 2160  
gtttttccat aggtccgccc cccctgacga gcatcacaaa aatcgacgct caagtcagag 2220  
gtggcgaaac ccgacaggac tataaagata ccaggcggtt cccctggaa gctccctcgt 2280  
gcgctctcct gttccgaccc tgccgcttac cggataacctg tcgcctttc tcccttcggg 2340  
aagcgtggcg ctttctcata gtcacgctg taggtatctc agttcgggtg aggtcgttcg 2400  
ctccaagctg ggctgtgtgc acgaaccccc cgttcagccc gaccgctgcg ccttatccgg 2460  
taactatcgt cttgagtcca acccgtaag acacgactta tcgccactgg cagcagccac 2520  
tggtaacagg attagcagag cgaggtatgt aggcgggtgct acagagttct tgaagtgggtg 2580  
gcctaactac ggctacacta gaaggacagt atttgggtatc tgcgctctgc tgaagccagt 2640  
taccttcgga aaaagagttg gtagctcttg atccggcaaa caaaccaccg ctggtagcgg 2700  
tgggtttttt gtttgcaagc agcagattac gcgcagaaaa aaaggatctc aagaagatcc 2760  
tttgatcttt tctacggggg ctgacgctca gtggaacgaa aactcacgtt aagggtttt 2820

```

ggatcatgaga ttatcaaaaa ggatcttcac ctagatcctt ttaaattaaa aatgaagttt 2880
taaatcaatc taaagtatat atgagtaaac ttgggtctgac agttaccaat gcttaatcag 2940
tgaggcacct atctcagcga tctgtctatt tcgttcaccc atagttgcct gactccccgt 3000
cgtgtagata actacgatac gggagggtt accatctggc cccagtgtg caatgatacc 3060
gcgagaccca cgctcaccgg ctccagattt atcagcaata aaccagccag ccggaagggc 3120
cgagcgaga agtggtcctg caactttatc cgcctccatc cagtctatta attgttgccg 3180
ggaagctaga gtaagtagtt cgccagttaa tagtttgccg aacgttggtg ccattgctac 3240
aggcatcgtg gtgtcacgct cgctgtttgg tatggcttca ttcagctccg gttcccaacg 3300
atcaaggcga gttacatgat ccccatgtt gtgcaaaaaa gcggttagct ccttcgggtcc 3360
tccgatcgtt gtcagaagta agttggccgc agtggtatca ctcatggtta tggcagcact 3420
gcataattct cttactgtca tgccatccgt aagatgcttt tctgtgactg gtgagtactc 3480
aaccaagtca ttctgagaat agtgatgctg gcgaccgagt tgctcttgcc cggcgtcaat 3540
acgggataat accgcgccac atagcagaac tttaaaagtg ctcatcattg gaaaacgttc 3600
ttcggggcga aaactctcaa ggatcttacc gctgttgaga tccagttcga tgtaaccac 3660
tcgtgcaccc aactgatctt cagcatcttt tactttcacc agcgtttctg ggtgagcaaa 3720
aacaggaagg caaatgccg caaaaaaggg aataaggcg acacggaaat gttgaatact 3780
catactcttc ctttttcaat attattgaag catttatcag ggttattgtc tcatgagcgg 3840
atacatattt gaatgtattt agaaaaataa acaaataggg gttccgcgca catttccccg 3900
aaaagtgcc cctgacgtct aagaaacat tattatcatg acattaacct ataaaaatag 3960
gcgtatcacg aggccctttc gtc 3983

```

&lt;210&gt; 2

&lt;211&gt; 16

&lt;212&gt; DNA

&lt;213&gt; Artificial Sequence

&lt;220&gt;

&lt;223&gt; Description of Artificial Sequence: primer MDB286

&lt;400&gt; 2

ngtcgaswga nawgaa

16

&lt;210&gt; 3

&lt;211&gt; 22

&lt;212&gt; DNA

&lt;213&gt; Artificial Sequence

&lt;220&gt;

&lt;223&gt; Description of Artificial Sequence: primer MDB391

&lt;400&gt; 3

tgatacaag catggtggat gg

22

&lt;210&gt; 4

&lt;211&gt; 21

&lt;212&gt; DNA

&lt;213&gt; Artificial Sequence

<220>

<223> Description of Artificial Sequence: primer MDB411

<400> 4

aggcatgccg ctgaaatcac c

21

<210> 5

<211> 21

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: primer MDB420

<400> 5

ggtttcgctc atgtgttgag c

21

<210> 6

<211> 1073

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: sequence  
comprising a 5' flanking region of GAT-ZM1

<220>

<221> misc\_feature

<222> (1)..(341)

<223> plant DNA

<220>

<221> misc\_feature

<222> (342)..(1073)

<223> T-DNA of pUC/Ac

<220>

<221> primer\_bind

<222> (286)..(307)

<223> binding site of primer COR17

<220>

<221> primer\_bind

<222> Complement((466)..(487))

<223> binding site of primer COR18

<400> 6

cgtcgagtga gatgaagtca cgacggggac tgactgcacc gtcgtctcag gtacgagggg 60

```

gacgtccagc aagcgtttcg cgagcvtgcc ggcgtcgtcc gtttgctcgg gattggcgtg 120
tcgcggggag acvgcvchcg tctttgtctc aaacvmgagg tcgatgcccg acgcgcccc 180
cgttggggcg ctggcgccgt cgactcgatc gacagccgac gaggcgctgc ctctgcttg 240
accttggttg cctgcctcc tcctccgtcg gcgggggaga ggacggggtg agctcgaatg 300
ttgttcttcc accacgcggg gaagacgtcg tcgattccac ctcatactc ttctttttc 360
aatattattg aagcatttat cagggttatt gtctcatgag cggatacata ttggaatgta 420
tttagaaaaa taaacaaata ggggttccgc gcacatttcc ccgaaaagtg ccacctgacg 480
tctaagaaac cattattatc atgacattaa cctataaaaa taggcgtatc acgaggccct 540
ttcgtctcgc gcgtttcggg gatgacgggtg aaaacctctg acacatgcag ctcccggaga 600
cggtcacagc ttgtctgtaa gcggatgccg ggagcagaca agcccgtcag ggcgcgtcag 660
cgggtggttg cgggtgtcgg ggctggctta actatgcggc atcagagcag attgtactga 720
gagtgcacca tatgcggtgt gaaataccgc acagatgcgt aaggagaaaa taccgcatca 780
ggcgccattc gccattcagg ctgcgcaact gttgggaagg gcgatcgggtg cgggcctctt 840
cgctattacg ccagctggcg aaagggggat gtgctgcaag gcgattaagt tgggtaacgc 900
cagggttttc ccagtcacga cgttgtaaaa cgacggccag tgccaagctt gaattcgagc 960
tcggtaccca ctggattttg gttttaggaa ttagaaattt tattgataga agtattttac 1020
aaatacaaat acatactaag ggtttcttat atgctcaaca catgagcgaa acc 1073

```

&lt;210&gt; 7

&lt;211&gt; 22

&lt;212&gt; DNA

&lt;213&gt; Artificial Sequence

&lt;220&gt;

&lt;223&gt; Description of Artificial Sequence: primer MDB439

&lt;400&gt; 7

ctcatggtta tggcagcact gc

22

&lt;210&gt; 8

&lt;211&gt; 23

&lt;212&gt; DNA

&lt;213&gt; Artificial Sequence

&lt;220&gt;

&lt;223&gt; Description of Artificial Sequence: primer VDS44

&lt;400&gt; 8

ctgtcatgcc atccgtaaga tgc

23

&lt;210&gt; 9

&lt;211&gt; 20

&lt;212&gt; DNA

&lt;213&gt; Artificial Sequence

&lt;220&gt;

&lt;223&gt; Description of Artificial Sequence: primer MDB522



<400> 9  
tgctttgaag acgtggttg

20

<210> 10  
<211> 484  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: sequence  
comprising a 3' flanking region of GAT-ZM1

<220>  
<221> misc\_feature  
<222> (1)..(342)  
<223> T-DNA of pUC/Ac

<220>  
<221> misc\_feature  
<222> (343)..(484)  
<223> plant DNA

<400> 10  
tgctttgaag acgtggttg aacgtcttct tttccacga tgctcctcgt ggggtgggggt 60  
ccatcttttg gaccactgtc ggcagaggca tcttcaacga tggcctttcc tttatcgcaa 120  
tgatggcatt tgtaggagcc accttccttt tctactatct tcataataaa gtgacagata 180  
gctgggcaat ggaatccgag gaggtttccg gatattaccc tttgttgaaa agtctcaatt 240  
gccctttggt cttctgagac tgtatctttg atatttttg agtagacaag cgtgtcgtgc 300  
tccacatgt tgacgaagat tttcttcttg tcattgagtc gttccgccat tgctcgtgc 360  
gcacggcggt ggaaggagta tcatgtcgta gctgccgtca agctccagat gggcagtctc 420  
cagcaacctc tccggcccgg gacggtgctc cgttcggga gtcttgagtt catctcactc 480  
gacc 484

<210> 11  
<211> 22  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: primer COR17

<400> 11  
gggtgagctc gaatgttggt ct

22

<210> 12  
<211> 22  
<212> DNA  
<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: primer COR18

<400> 12

tcttagacgt caggtggcac tt

22

<210> 13

<211> 22

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: primer COR15

<400> 13

agcgtcaagg atcattggtg tc

22

<210> 14

<211> 22

<212> DNA

<213> Artificial Sequence

<220>

<223> Description of Artificial Sequence: primer COR16

<400> 14

ggccaagttc agcataagct gt

22

(12) INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(19) World Intellectual Property Organization  
International Bureau



(43) International Publication Date  
19 July 2001 (19.07.2001)

PCT

(10) International Publication Number  
**WO 01/51654 A3**

(51) International Patent Classification<sup>7</sup>: C12Q 1/68.  
C12N 15/82

(21) International Application Number: PCT/EP01/00306

(22) International Filing Date: 9 January 2001 (09.01.2001)

(25) Filing Language: English

(26) Publication Language: English

(30) Priority Data:  
09/481,049 11 January 2000 (11.01.2000) US

(71) Applicant: AVENTIS CROPSOURCE N.V. [BE/BE];  
Jozef Plateastraat 22, B-9000 Gent (BE).

(72) Inventor: DE BEUCKELEER, Marc; Hekers 36,  
B-9052 Zwijnaarde (BE).

(81) Designated States (national): AE, AG, AL, AM, AT, AU,  
AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CR, CU, CZ,

DE, DK, DM, DZ, EE, ES, FI, GB, GD, GE, GH, GM, HR,  
HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR,  
LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ,  
NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM,  
TR, TT, TZ, UA, UG, UZ, VN, YU, ZA, ZW.

(84) Designated States (regional): ARIPO patent (GH, GM,  
KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZW), Eurasian  
patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European  
patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE,  
IT, LU, MC, NL, PT, SE, TR), OAPI patent (BF, BJ, CF,  
CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG).

**Published:**

— with international search report

(88) Date of publication of the international search report:  
31 January 2002

*For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.*



**WO 01/51654 A3**

(54) Title: METHODS AND KITS FOR IDENTIFYING ELITE EVENT GAT-ZM1 IN BIOLOGICAL SAMPLES

(57) Abstract: Tools are provided which allow rapid and unequivocal identification elite event GAT-ZM1 in biological samples.

# INTERNATIONAL SEARCH REPORT

International Application No  
PC1/EP 01/00306

**A. CLASSIFICATION OF SUBJECT MATTER**  
IPC 7 C12Q1/68 C12N15/82

According to International Patent Classification (IPC) or to both national classification and IPC

**B. FIELDS SEARCHED**

Minimum documentation searched (classification system followed by classification symbols)  
IPC 7 C12Q C12N

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

EPO-Internal, SEQUENCE SEARCH, WPI Data, PAJ, MEDLINE, BIOSIS, EMBASE, CHEM ABS Data

**C. DOCUMENTS CONSIDERED TO BE RELEVANT**

Category	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No
X	WO 95 22625 A (AFFYMAX TECH NV ; CRAMER ANDREAS (US); STEMMER WILLEM P C (US)) 24 August 1995 (1995-08-24)	12, 13, 15, 22-25
Y	* see especially page 84, line 36 (primer 2650) * abstract	1-8, 10, 11, 16-21, 26, 27
Y	US 5 489 520 A (CHAMBERS SHERYL A ET AL) 6 February 1996 (1996-02-06)  the whole document	1-8, 10, 11, 16-21, 26, 27
A	WO 98 30701 A (MEYER TERRY; JAYNE SUSAN ; BARBOUR ERIC (US)) 16 July 1998 (1998-07-16) the whole document	

-/--

☒ Further documents are listed in the continuation of box C

☒ Patent family members are listed in annex

\* Special categories of cited documents:

- \*A\* document defining the general state of the art which is not considered to be of particular relevance
- \*E\* earlier document but published on or after the international filing date
- \*L\* document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
- \*O\* document referring to an oral disclosure, use, exhibition or other means
- \*P\* document published prior to the international filing date but later than the priority date claimed

- \*T\* later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
- \*X\* document of particular relevance: the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
- \*Y\* document of particular relevance: the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.
- \*G\* document member of the same patent family

Date of the actual completion of the international search

12 October 2001

Date of mailing of the international search report

26/10/2001

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2  
NL - 2280 HV Rijswijk  
Tel. (+31-70) 340-2040, Tx. 31 651 epo nl,  
Fax: (+31-70) 340-3016

Authorized officer

Knehr, M

## INTERNATIONAL SEARCH REPORT

International Application No

PC1/EP 01/00306

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT		
Category	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No
A	EP 0 275 957 A (HOECHST AG) 27 July 1988 (1988-07-27) cited in the application the whole document ---	
A	LIU Y-G ET AL.: "Efficient isolation and mapping of Arabidopsis thaliana T-DNA insert junctions by thermal asymmetric interlaced PCR" THE PLANT JOURNAL, vol. 8, no. 3, 1995, pages 457-463, XP002180068 cited in the application the whole document ---	
A	THOMPSON C J ET AL.: "Characterization of the herbicide-resistance gene bar from Streptomyces hygrosopicus" THE EMBO JOURNAL, vol. 6, no. 9, 1987, pages 2519-2523, XP002180069 the whole document ---	
T	WO 01 31042 A (AVENTIS CROPS SCIENCE N V) 3 May 2001 (2001-05-03) * see especially claims 12-20 * the whole document -----	1-27

## INTERNATIONAL SEARCH REPORT

Information on patent family members

International Application No.

PCT/EP 01/00306

Patent document cited in search report		Publication date	Patent family member(s)	Publication date
WO 9522625	A	24-08-1995	US 5605793 A	25-02-1997
			AU 703264 B2	25-03-1999
			AU 2971495 A	04-09-1995
			CA 2182393 A1	24-08-1995
			CN 1145641 A	19-03-1997
			DE 752008 T1	05-04-2001
			DE 934999 T1	05-04-2001
			EP 1094108 A2	25-04-2001
			EP 0752008 A1	08-01-1997
			EP 0934999 A1	11-08-1999
			JP 10500561 T	20-01-1998
			JP 2001057893 A	06-03-2001
			WO 9522625 A1	24-08-1995
			US 6180406 B1	30-01-2001
			US 6132970 A	17-10-2000
			US 6287861 B1	11-09-2001
			US 6291242 B1	18-09-2001
			US 6277638 B1	21-08-2001
			US 6297053 B1	02-10-2001
			US 5830721 A	03-11-1998
			US 5811238 A	22-09-1998
			US 6117679 A	12-09-2000
			US 5837458 A	17-11-1998
US 5489520	A	06-02-1996	US 5550318 A	27-08-1996
			US 5919675 A	06-07-1999
			US 5969213 A	19-10-1999
			US 5874265 A	23-02-1999
			AT 160174 T	15-11-1997
			AU 644097 B2	02-12-1993
			AU 6290390 A	11-03-1991
			BR 9007588 A	30-06-1992
			CA 2064761 A1	10-02-1991
			DE 69031711 D1	18-12-1997
			DE 69031711 T2	09-04-1998
			EP 0485506 A1	20-05-1992
			EP 0814166 A2	29-12-1997
			ES 2110417 T3	16-02-1998
			HU 60781 A2	28-10-1992
			JP 5501352 T	18-03-1993
			WO 9102071 A2	21-02-1991
WO 9830701	A	16-07-1998	AU 729122 B2	25-01-2001
			AU 5734198 A	03-08-1998
			EP 0954584 A1	10-11-1999
			US 6096947 A	01-08-2000
			WO 9830701 A1	16-07-1998
			ZA 9800260 A	13-07-1998
EP 0275957	A	27-07-1988	DE 3737918 A1	02-03-1989
			AT 86300 T	15-03-1993
			AU 609082 B2	26-04-1991
			AU 1061988 A	28-07-1988
			CA 1321364 A1	17-08-1993
			CN 87100603 A	10-08-1988
			CN 88100322 A .B	10-08-1988
			DE 3878699 D1	08-04-1993
			DK 23988 A	22-07-1988

Form PCT/ISA/210 (patent family annex) (July 1992)

# INTERNATIONAL SEARCH REPORT

Information on patent family members

International Application No

PCT/EP 01/00306

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
EP 0275957	A	EP 0275957 A2	27-07-1988
		ES 2054708 T3	16-08-1994
		FI 880216 A	22-07-1988
		GR 3007859 T3	31-08-1993
		HU 215079 B	28-09-1998
		JP 2993964 B2	27-12-1999
		JP 63273479 A	10-11-1988
		JP 3062125 B2	10-07-2000
		JP 10080289 A	31-03-1998
		JP 3093686 B2	03-10-2000
		JP 10080278 A	31-03-1998
		NZ 223227 A	28-06-1989
		US 5276268 A	04-01-1994
		ZA 8800368 A	05-07-1988
WO 0131042	A	03-05-2001	
		AU 1390001 A	08-05-2001
		WO 0131042 A2	03-05-2001

